



Elastic Bumper Kit

ME-8998

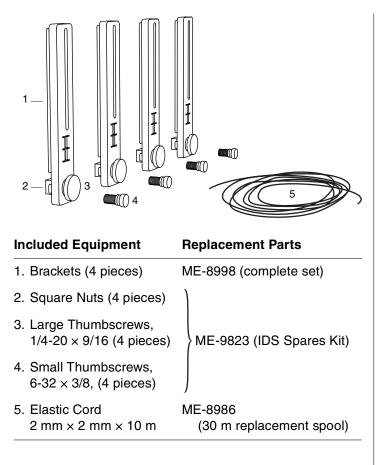

Introduction

The Elastic Bumper Kit includes parts that can be used as a pair of elastic bumpers or for mounting four photogates on a PASCO dynamics track.

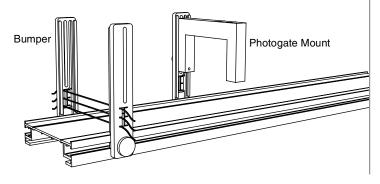


To Mount Bracket on Track

The brackets attach to the T-slots on the sides of a dynamics track.

- 1. Loosen (but do not remove) the square nut on the thumb-screw.
- 2. Slide the nut into the T-slot at the end of the track.
- **3.** Slide the bracket along the track to the desired location.
- **4.** Tighten the thumbscrew to secure the bracket.

For Use as Bumpers


Using the four included brackets and some of the included elastic cord, you can set up two bumpers on a track. This type of bumper is ideal for use with a motion sensor because the sensor can "see through" the thin elastic cord, and the

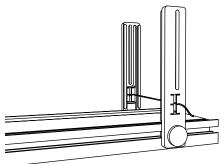

bumper protects the sensor from being hit by a cart.


- 1. Mount two brackets on either side of the track as illustrated.
- 2. Cut off a 15 cm length of elastic cord.
- **3.** Stretch the cord through the lower slots in the brackets. Secure the cord by pulling it into the notches.

- 4. For heavy carts or fast collisions, add one or two more pieces of elastic cord.
- 5. If you are using a motion sensor, place it within 15 cm of the bumper. This will allow the sensor to "see" through the bumper.

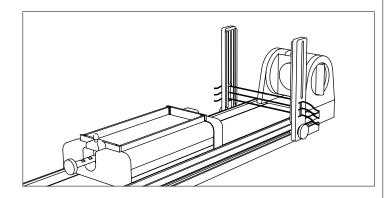

For Use as Photogate Mount


- 1. Mount a bracket on the track at the desired location.
- 2. Remove the rod-mounting clamp from the photogate.

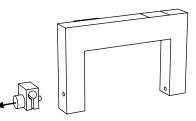

3. Insert a small thumbscrew through the upper slot of the bracket and into the mounting hole in the photogate.

Note: Use the shorter thumbscrew (3/8 inch, 10 mm) included with the bumper kit, not the longer screw that came with the photogate.

- 4. Adjust the height of the photogate.
- 5. Tighten the thumbscrew to secure the photogate.

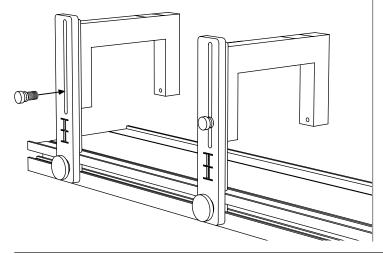

Technical Support

For assistance with any PASCO product, contact PASCO at:

Address:	PASCO scientific
	10101 Foothills Blvd.
	Roseville, CA 95747-7100
Phone:	916-786-3800 (worldwide)
	800-772-8700 (U.S.)
Fax:	(916) 786-7565
Web:	www.pasco.com
Email:	support@pasco.com

Limited Warranty For a description of the product warranty, see the PASCO catalog.

Copyright The PASCO scientific 012-10279A *Elastic Bumper Kit Instruction Sheet* is copyrighted with all rights reserved. Permission is granted to non-profit educational institutions for reproduction of any part of this manual, providing the reproductions are used only in their laboratories and classrooms, and are not sold for profit. Reproduction under any other circumstances, without the written consent of PASCO scientific, is prohibited.

Trademarks PASCO and PASCO scientific are trademarks or registered trademarks of PASCO scientific, in the United States and/or in other countries. All other brands, products, or service names are or may be trademarks or service marks of, and are used to identify, products or services of, their respective owners. For more information visit www.pasco.com/legal.